



Microscope objectives

Our objectives help you focus on yours

Nikon is a leader in the development and manufacture of optical and digital imaging technology for advanced science and clinical research. With over a 100-year history of optical excellence, Nikon is committed to accelerating innovation in science and clinical imaging to improve healthcare and provide a better quality of life.



Nikon's first microscope, released in 1925

The switch from traditional film photography to digital imaging was a major milestone in the field of microscopy, opening up new possibilities in both application and technology. Introduction of digital imaging spurred significant technological changes including the development of objectives with enhanced optical quality and functionality to meet the new demands. Objective lenses are arguably the most important element in the microscope and Nikon continues to invest heavily in the development of objectives to meet the changing demands of science. Explore some of Nikon's newest developments in high-performance objectives in this brochure.



A vertical beam of light illuminates a dark industrial setting, likely a glass melting process. The beam is bright and focused, creating a strong contrast with the surrounding dark environment. The background shows structural elements of a factory or laboratory, including beams and pipes, all dimly lit.


In tireless pursuit of the highest quality

Each Nikon microscope objective is precision-crafted to provide the highest level of clarity and overall optical performance. World-class Nikon objectives, including the renowned CFI60 infinity optical system, deliver brilliant images of breathtaking sharpness and clarity, from ultralow to the highest magnifications.

The initial melting process after blending glass raw materials

It Starts with the Glass

Nikon has been developing optical glass since its inception in 1917, and to this day, wholly owns and formulates all of its glass. Optical glass starts as an ingot (shown on right) which is formed by blending rare earth elements and repeated melting, shaping and slow cooling to achieve a target refractive index. The glass ingots are precision-cut, polished and coated to produce lens elements for the objective.

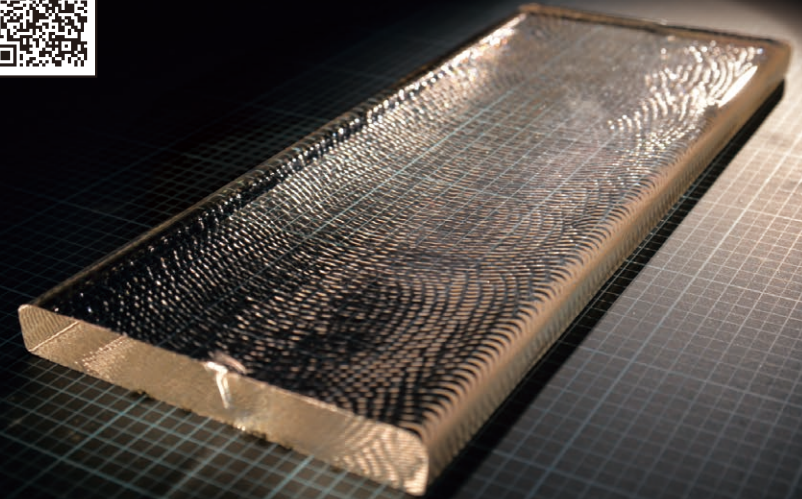


Mastering Excellence

The front lens of high-performance objectives is hand-polished by Nikon's most highly skilled experts (shown on left), a technique requiring more than a decade to master.

By controlling the entire manufacturing process from glass formulation to assembly and alignment of lens elements, Nikon ensures the highest quality and performance of its objectives.

Objectives Inspired by Your Science

Produced to exacting standards, Nikon's objectives provide exceptional detail and clarity. The highest level of image quality can be achieved whether it be for routine tasks or cutting-edge research.

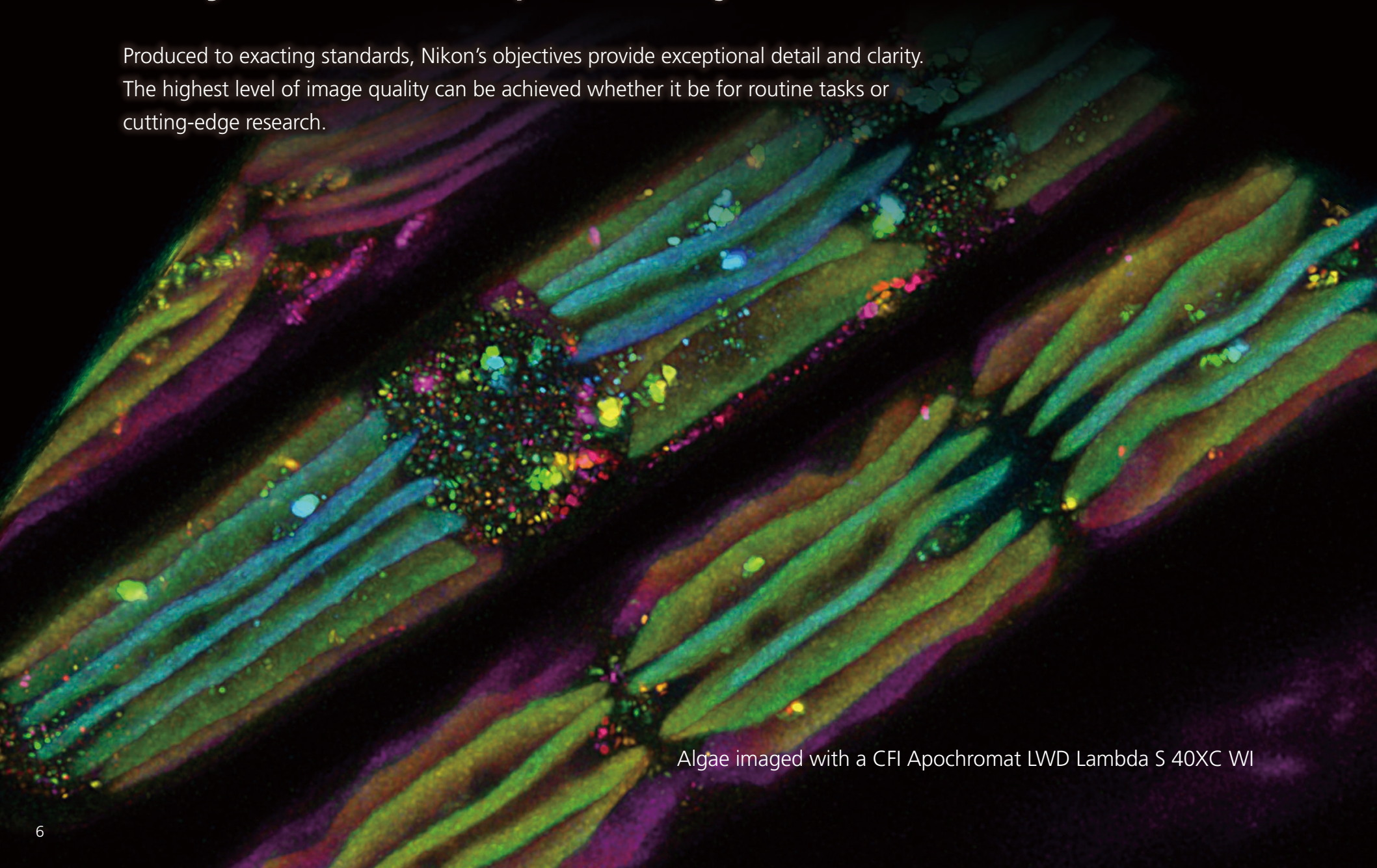

Algae imaged with a CFI Apochromat LWD Lambda S 40XC WI

Cleared Tissue Imaging




CFI90 20XC Glyc

Designed for deep imaging of cleared tissues

The unique CFI90 design results in an incredibly high N.A. while maintaining a large field of view and ultra-long working distance. Incorporates a correction collar for compensating for different refractive indices of clearing agents.

- NA: 1.00, WD: 8.20 mm
- Chromatic aberration correction: from visible to IR
- High-transmittance Nano Crystal Coat
- Correction collar for spherical aberration

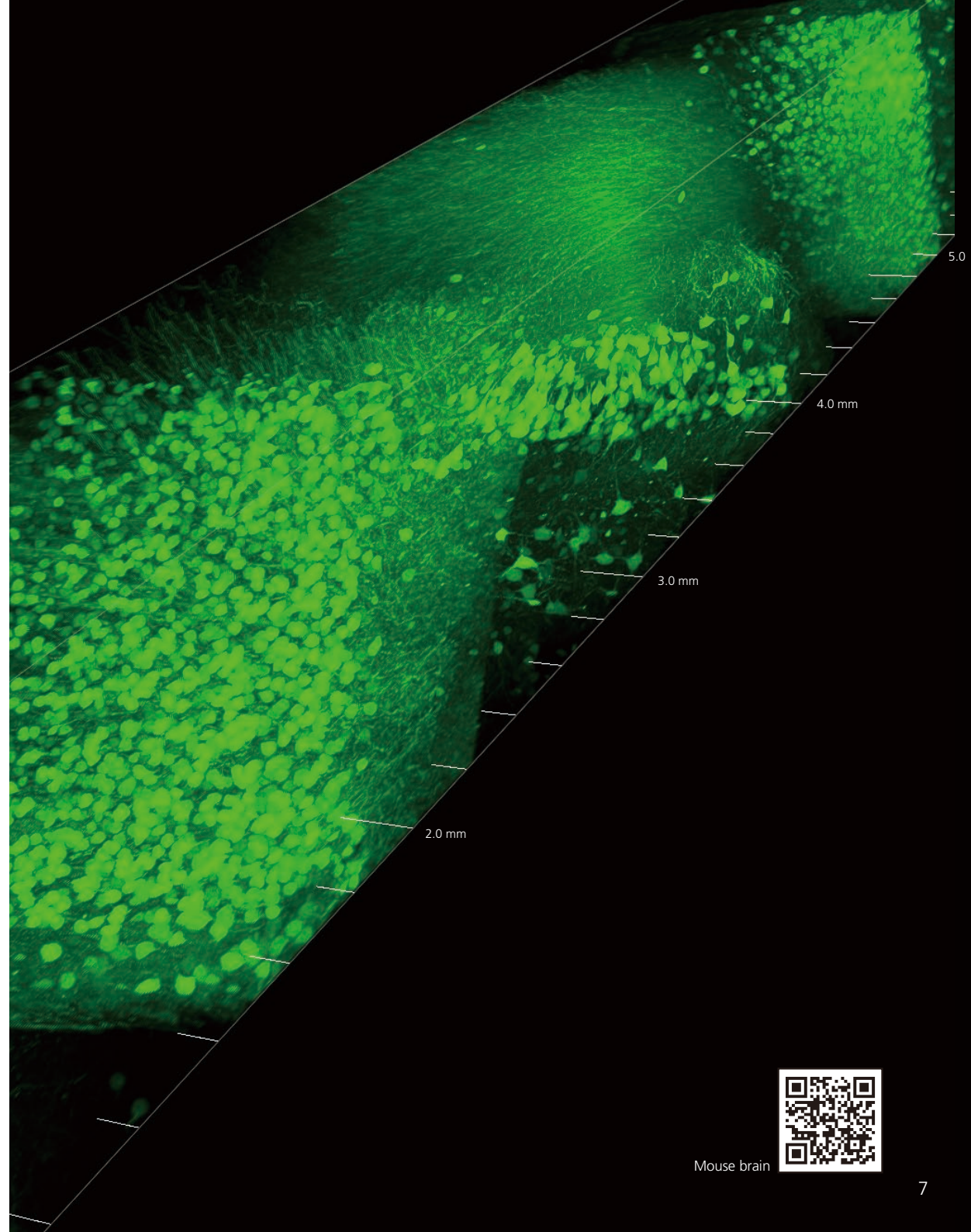

Mouse brain



Cleared Tissue Imaging




CFI Plan Apochromat 10XC Glyc

Compatible with a wide range of immersion media and clearing agents

In addition to water and immersion oil, this objective lens is compatible with a variety of tissue clearing agents. The lens also features chromatic aberration correction over a broad spectral range and is compatible with the Ti2 inverted microscope.

- NA: 0.50, WD: 5.50 mm (upright) / 2.00 mm (inverted)
- Chromatic aberration correction: from visible to IR
- High-transmittance Nano Crystal Coat
- Correction collar for spherical aberration

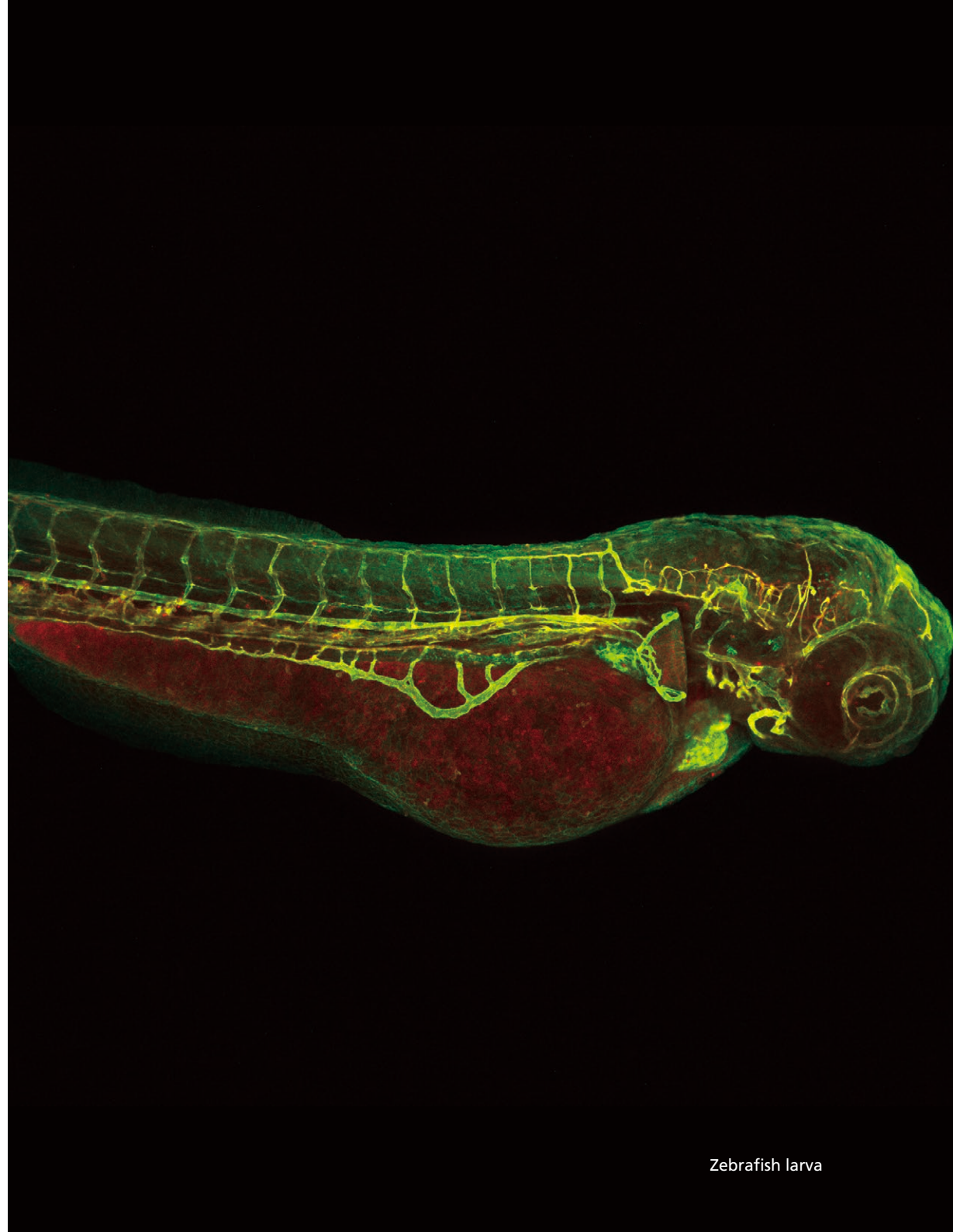

Zebrafish larva

Multiphoton Imaging




CFI75 Achromat 25XC W 1300

Developed for deep brain applications

Best in class objective for multiphoton and electrophysiology applications. Offers ultra-long working distance and incredibly high N.A. for exceptional results every time.

- NA: 1.10, WD: 2.00 mm
- Chromatic aberration correction: from visible to near IR
- High-transmittance Nano Crystal Coat
- Correction collar for spherical aberration

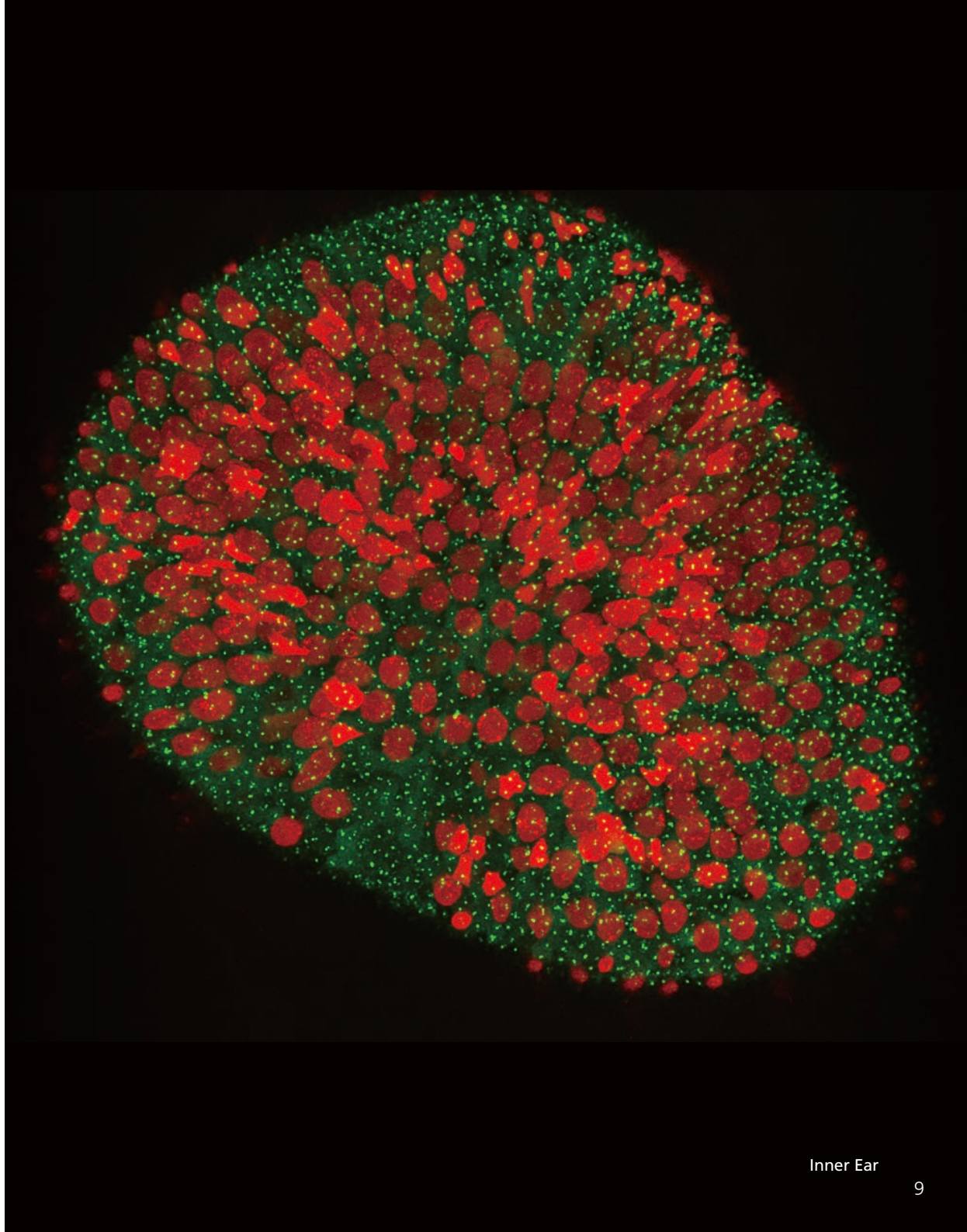

Confocal Imaging




CFI Apochromat LWD Lambda S 20XC WI

A versatile objective for a wide range of applications

Combines high N.A., long working distance and large field of view for imaging thick, live samples. Chromatic aberration correction from visible to IR for multiphoton imaging as well.

- NA: 0.95, WD: 0.95 mm
- Chromatic aberration correction: from visible to IR
- High-transmittance Nano Crystal Coat
- Correction collar for spherical aberration

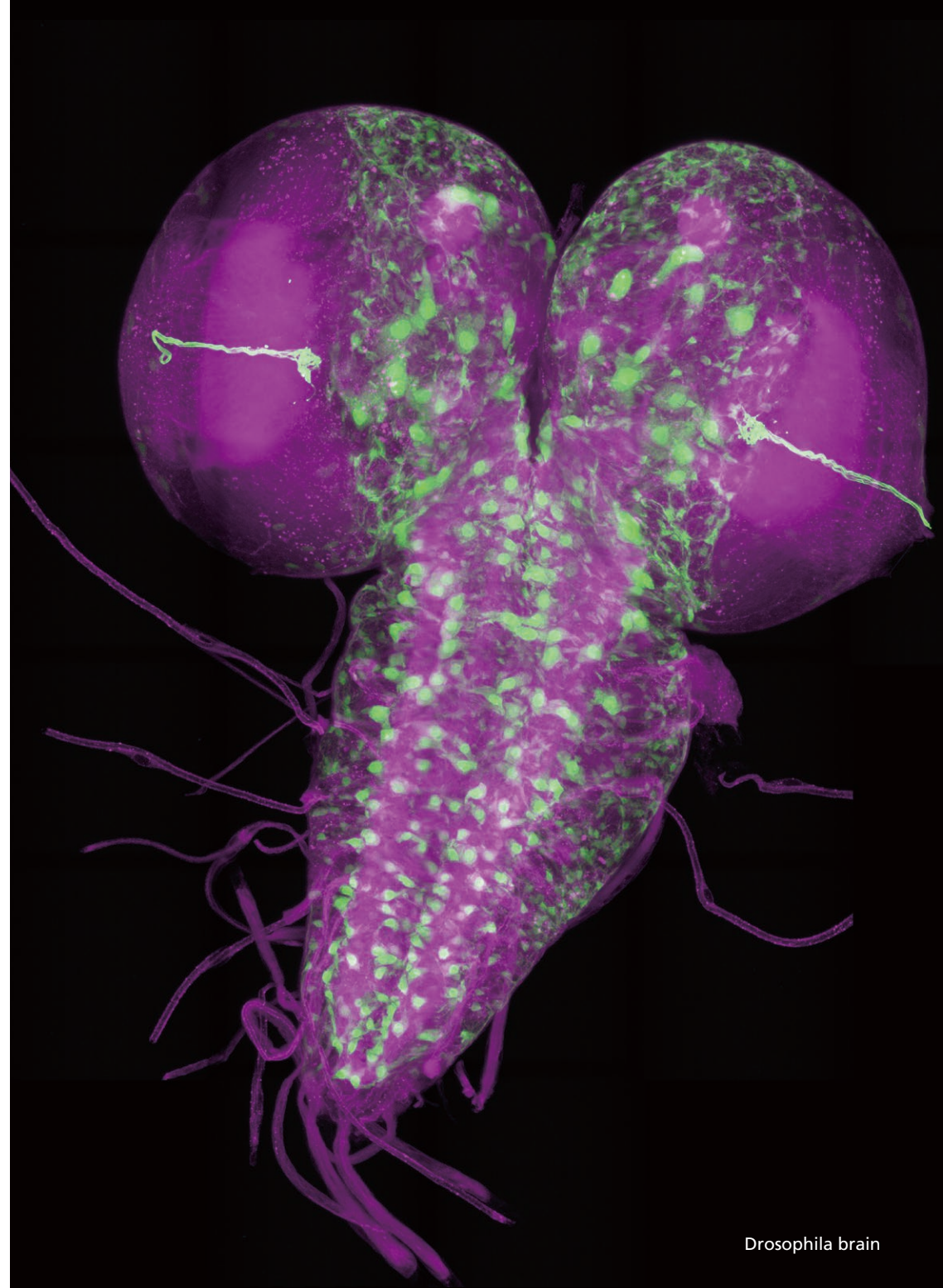

Drosophila brain

Confocal Imaging




CFI Apochromat Lambda S 40XC WI

Image live cells with stunning clarity

Boasting the highest N.A. of all 40X water immersion objectives, this lens provides incredible detail and clarity for multi-color live cell imaging applications. Can be used with an external phase contrast system (Ti2-E) for combining phase contrast and fluorescence imaging without attenuating transmission.

- NA: 1.25, WD: 0.18 mm
- Chromatic aberration correction: from visible to near IR
- High-transmittance Nano Crystal Coat
- Correction collar for spherical aberration

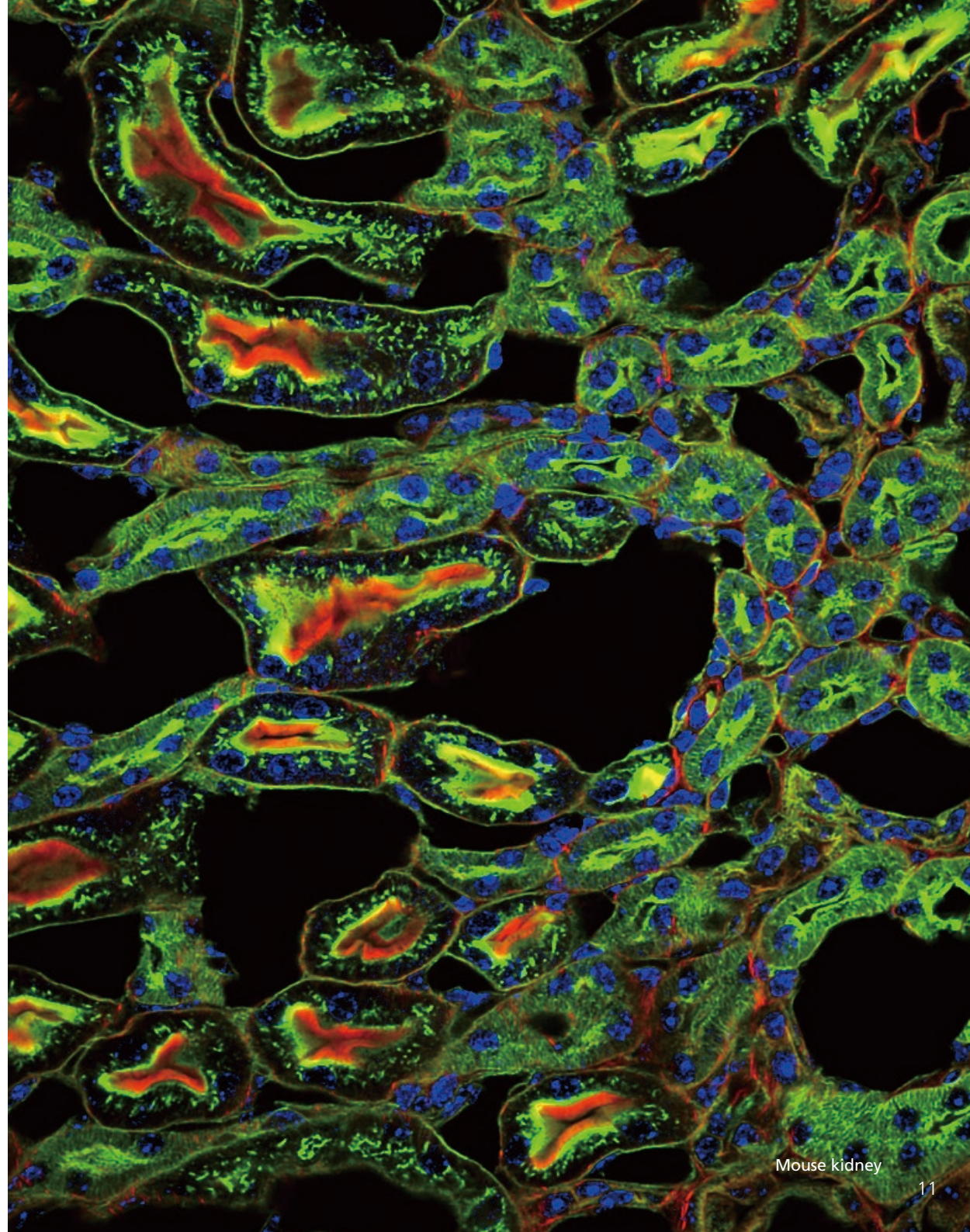

Confocal Imaging




CFI APOCHROMAT LWD Lambda S 40XC WI

A multi-purpose objective with long WD and high NA

This water immersion objective for inverted microscopes provides a long WD, high NA and high transmittance, and is effective for deep imaging with multiphoton microscopes. It is also compatible with DIC observation and external phase-contrast observation using the Ti2-E inverted microscope.

- NA: 1.15, WD: 0.60 mm
- Chromatic aberration correction: from visible to near IR
- High-transmittance Nano Crystal Coat
- Correction collar for spherical aberration

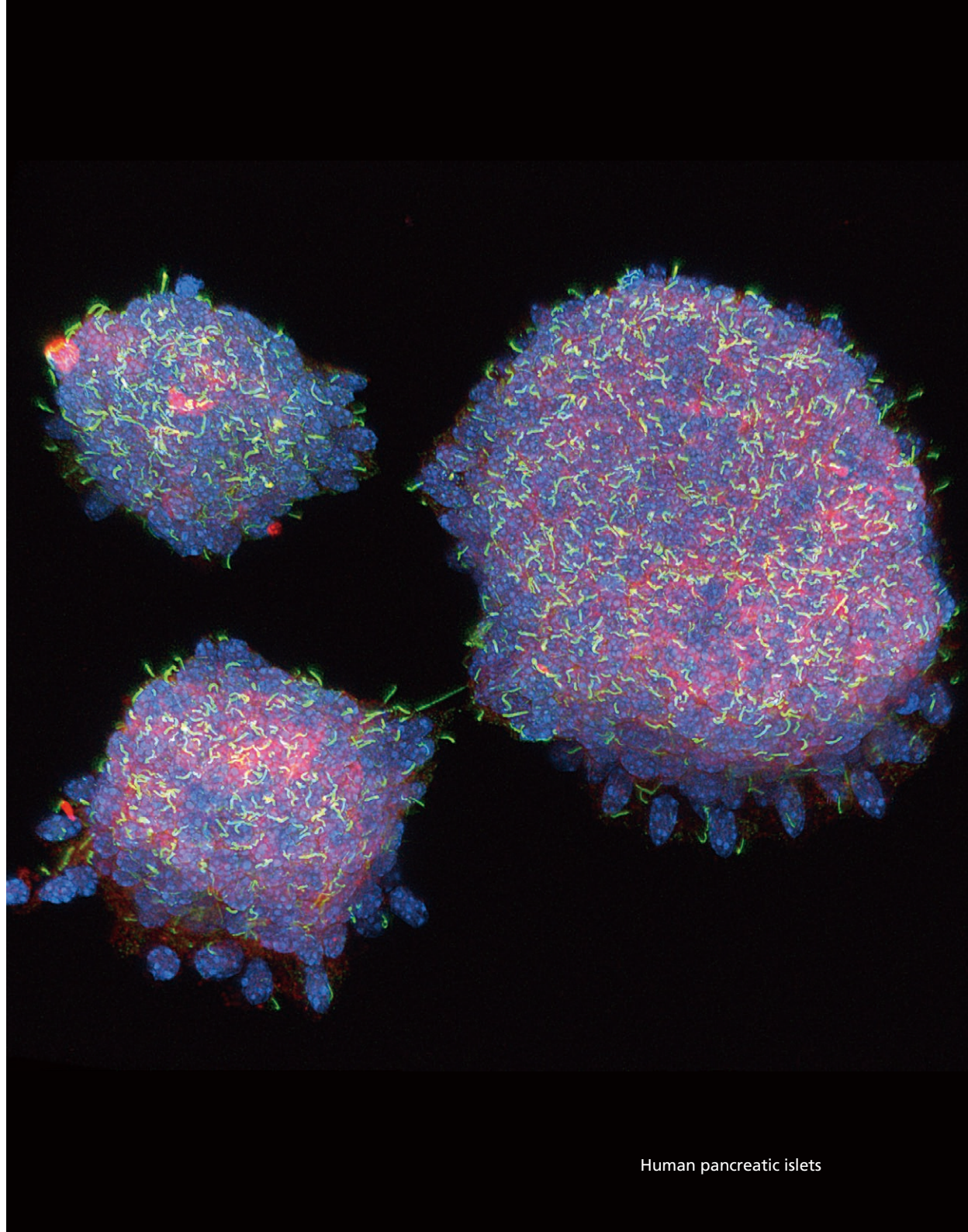

Human pancreatic islets

Confocal Imaging




CFI Plan Achromat IR 60XC WI

Incredible resolution and field flatness

This objective provides an NA of 1.27, the highest among 60X water immersion objectives, and achieves incredible field flatness. Corrects for chromatic aberration over a wide wavelength range up to IR, and supports various applications including multicolor live-cell confocal imaging, multiphoton imaging and laser tweezer applications.

- NA: 1.27, WD: 0.17 mm
- Chromatic aberration correction: from visible to IR
- High-transmittance Nano Crystal Coat
- Correction collar for spherical aberration

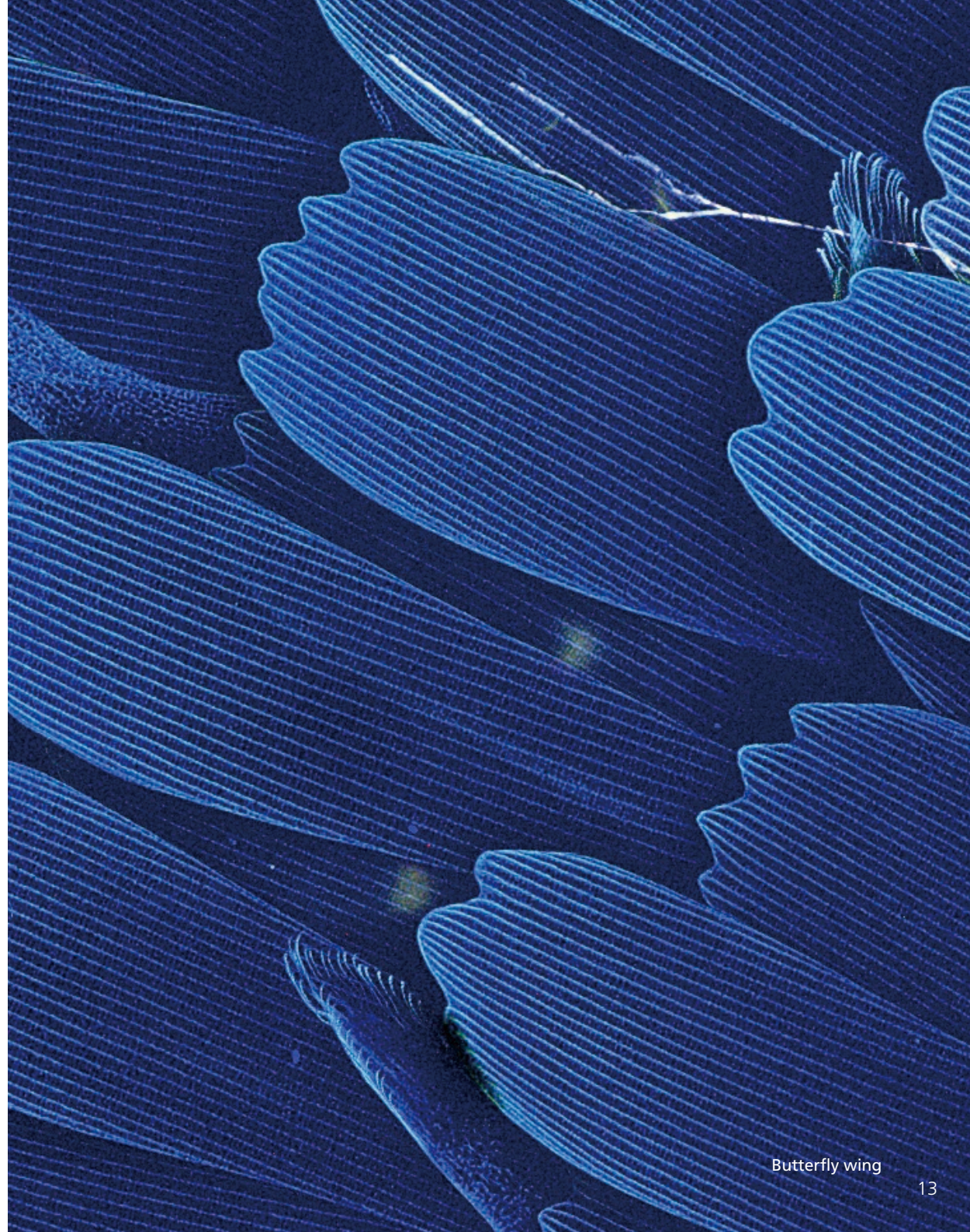

Super-Resolution Imaging




CFI SR HP Plan Apochromat Lambda S 100XC Sil

High-resolution silicone oil immersion objective for imaging live samples

Silicone oil more closely matches the refractive index of live cells compared to water or oil, thereby minimizing spherical aberration issues common to live cell imaging. In addition, the 100XC Sil lens achieves exceptional resolving power even at greater depths, making it well-suited for super-resolution imaging of thicker specimen.

- NA: 1.35, WD: 0.30 mm
- Chromatic aberration correction: visible range
- High-transmittance Nano Crystal Coat
- Correction collar for spherical aberration

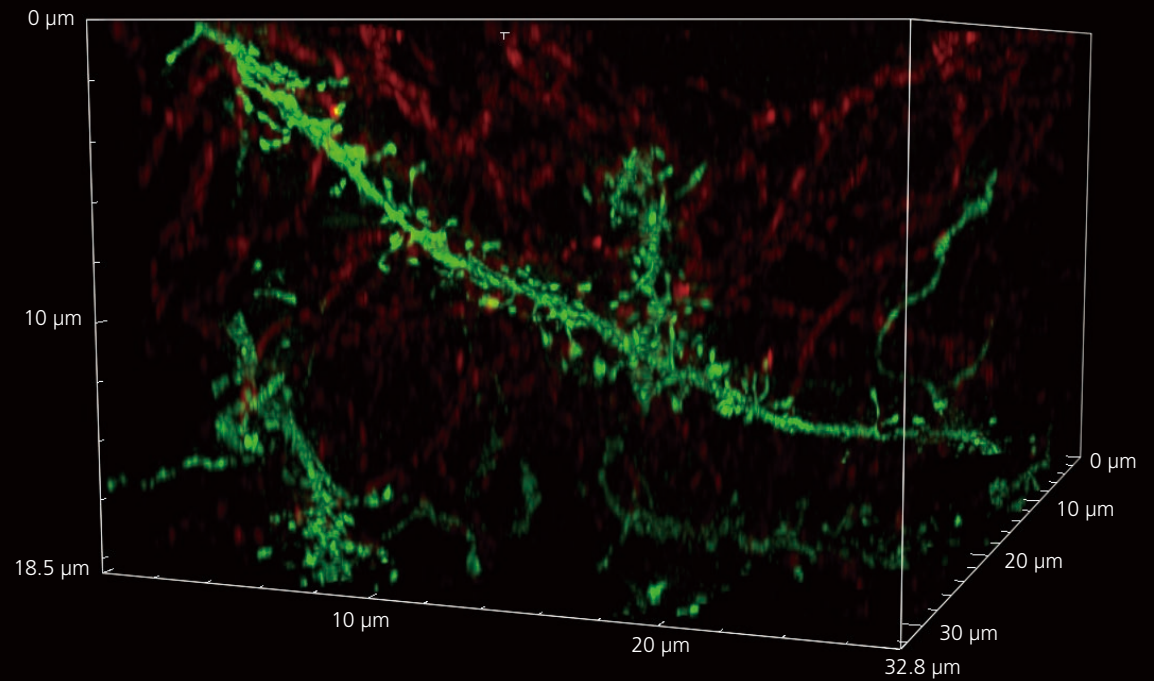

Dendrites


Super-Resolution Imaging



Auto Correction Collar

Quick and accurate spherical aberration correction

Achieving the highest quality point spread function is critical for super-resolution imaging. The ACC quickly moves lens elements in the objective to an optimal position based on the acquired point spread function to minimize spherical aberration.

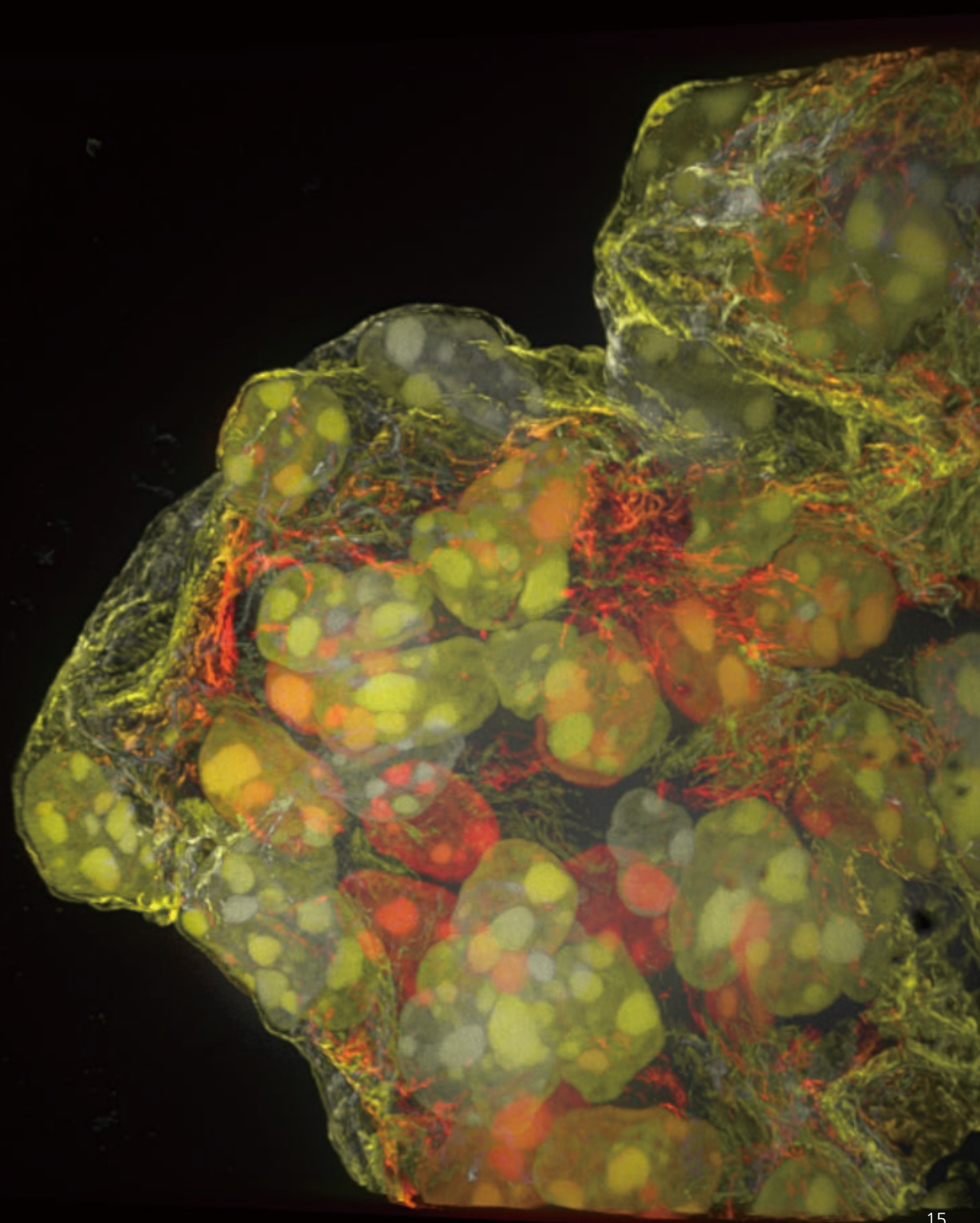

High-Content Imaging




CFI S Plan Fluor LWD ADM 20XC

High NA and long WD objective designed
for HC applications

This objective is compatible with thick plastic-bottom dishes and well-plates, and enables high-resolution phase-contrast and fluorescence observations. Its large FOV improves throughput of high-content applications.

- NA: 0.70, WD: 1.30-2.30 mm
- Chromatic aberration correction: visible range
- Correction collar for spherical aberration

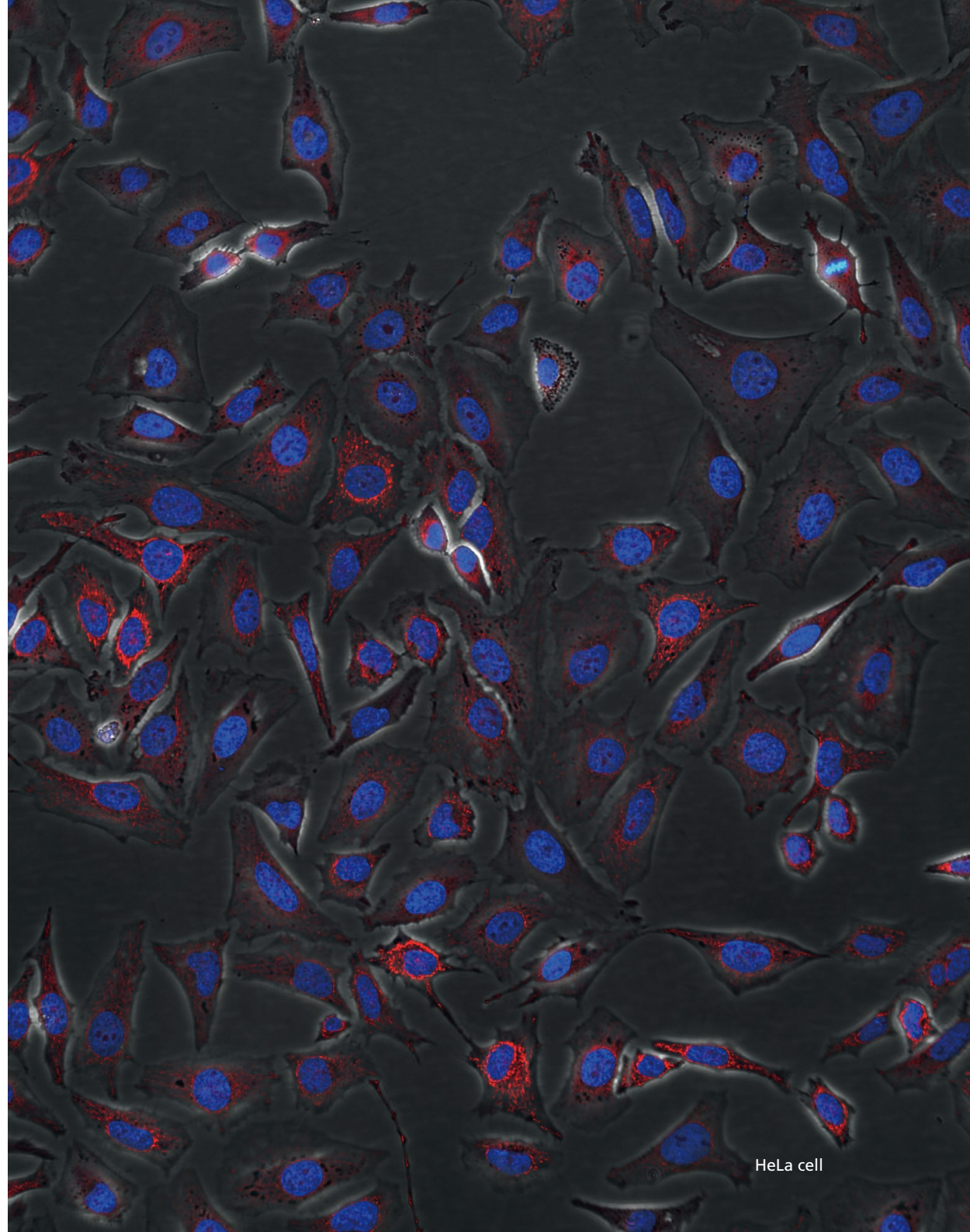

HeLa cell

Pathological Examination




CFI Plan Achromat Lambda 40XC

High-clarity, large field-of-view imaging for Clinical applications

High N.A. combined with Nano Crystal Coat results in bright, high-contrast images. Provides chromatic aberration correction over a large spectral range and superior image flatness from edge-to-edge over its 25 mm FOV. An optimal lens for both fluorescence and brightfield imaging.

- NA: 0.95, WD: 0.25-0.16 mm
- Chromatic aberration correction: from visible to near IR
- High-transmittance Nano Crystal Coat
- Correction collar for spherical aberration

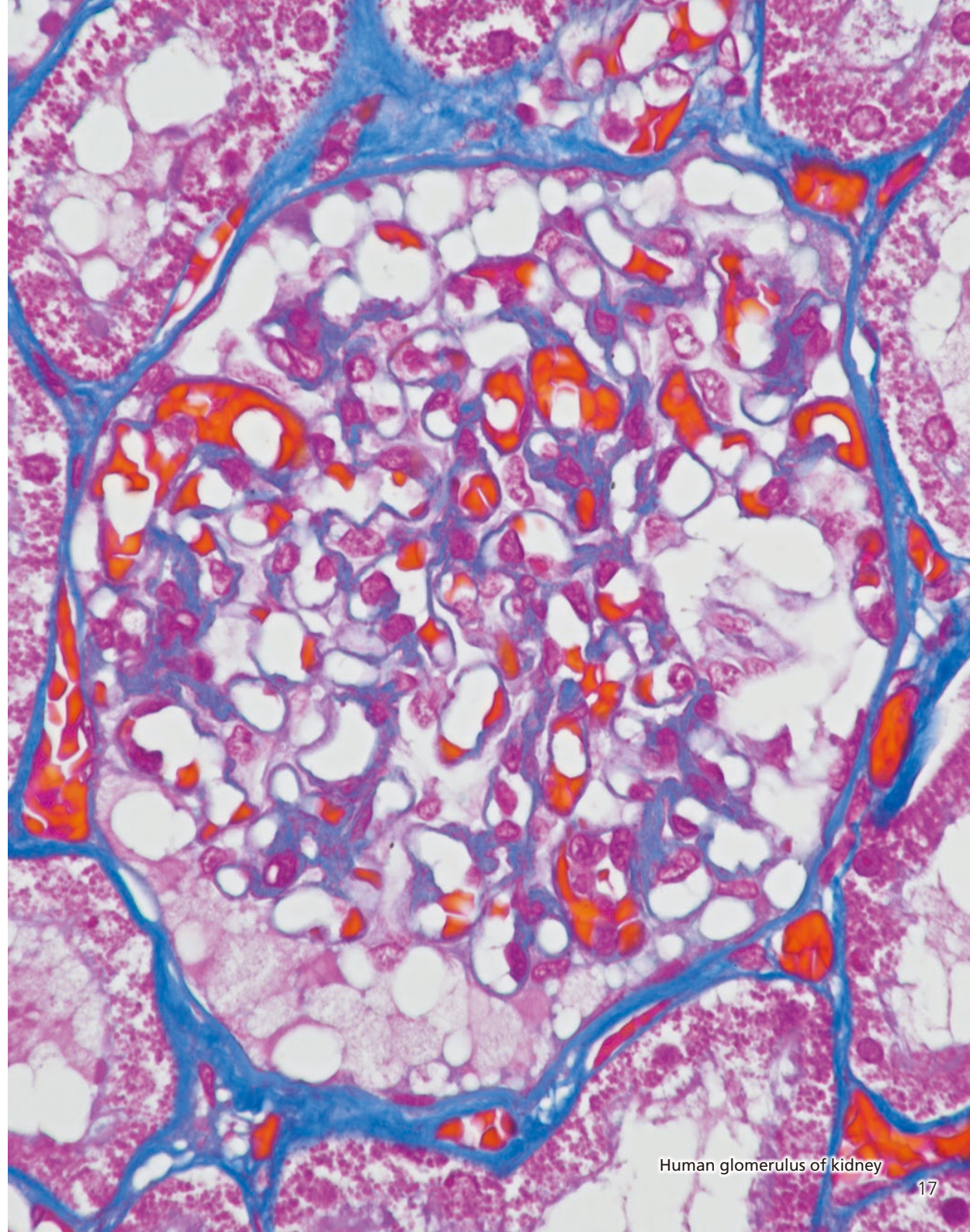


Human glomerulus of kidney

A tool for finding the right objective

Nikon's online Objective Selector tool enables you to quickly and easily find the right objective for your application. Refine your search based on application, technique, objective class, immersion type, etc. Specifications for multiple objectives can be displayed in a single window for easy comparison.

<https://www.microscope.healthcare.nikon.com/selectors/objectives>

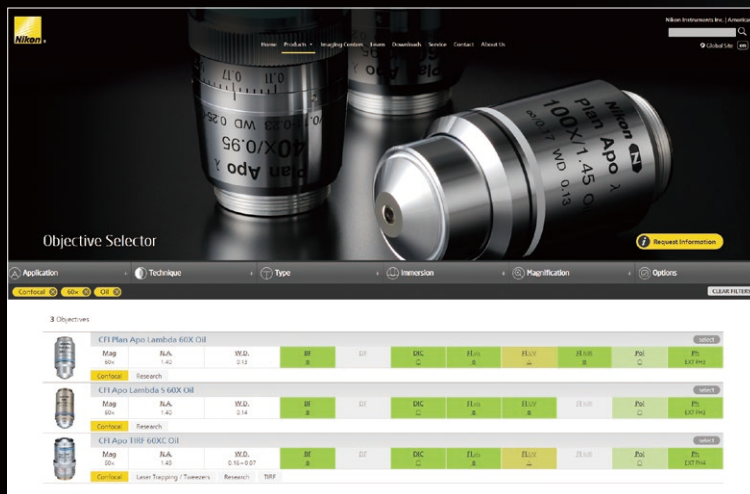

OBJECTIVE SELECTOR

As individual as you

Nikon offers a broad portfolio of objectives to meet your individual application needs. Explore some of the different objective series in the following pages.



CFI Plan Apochromat Lambda series

This series of objectives incorporates "Nano crystal coat" technology which dramatically improves light transmission over a wide spectral range, including near-IR. Chromatic aberration correction (up to near-IR), resolution and image flatness of these objectives are unrivaled, ensuring the highest quality images for a variety of applications.



CFI Plan Apochromat Lambda 2X, 4X, 10X, 20X, 40XC, 60XC, 60X Oil and 100X Oil



CFI Plan Apochromat DM Lambda 20X, 40XC, 60XC, 60X Oil and 100X Oil

Use	Model	NA	W.D. (mm)	Cover glass thickness	Brightfield	Darkfield	DIC	Phase contrast	Polarizing	Fluorescence			Ti2-E PFS
										Visible light	UV	NIR	
Brightfield (CFI Plan Apo)	Lambda 2X	0.10	8.50	—	⊙				○	⊙	△	⊙	
	Lambda 4X	0.20	20.00	—	⊙				○	⊙	△	⊙	●
	Lambda 10X	0.45	4.00	0.17	⊙	△	○		○	⊙	△	⊙	●
	Lambda 20X	0.75	1.00	0.17	⊙	○●	○		○	⊙	△	⊙	●
	Lambda 40XC	0.95	0.25-0.16	0.11-0.23	⊙	●	○		○	⊙	△	⊙	●
	Lambda 60XC	0.95	0.21-0.11	0.11-0.23	⊙	●	○		○	⊙	△	⊙	●
	Lambda 60X Oil	1.40	0.13	0.17	⊙			EXT PH3-60X	○	⊙	△	⊙	●
	Lambda 100X Oil	1.45	0.13	0.17	⊙			EXT PH3-100X EXT PH4-100X	○	⊙	△	⊙	●
	IR 60XC WI	1.27	0.18-0.16	0.15-0.19	⊙			EXT PH3-60x	○	○	△	⊙	●
NCG 100X Oil	1.40	0.16	0	⊙				○	⊙	△			
Phase contrast (CFI Plan Apo)	DM Lambda 20X	0.75	1.00	0.17	○	○●		⊙ PH2		○	△	○	●
	DM Lambda 40XC	0.95	0.25-0.16	0.11-0.23	○	●		⊙ PH2		○	△	○	●
	DM Lambda 60XC	0.95	0.21-0.11	0.11-0.23	○	●		⊙ PH2		○	△	○	●
	DM Lambda 60X Oil	1.40	0.13	0.17	○	●		⊙ PH3		○	△	○	●
	DM Lambda 100X Oil	1.45	0.13	0.17	○			⊙ PH3		○	△	○	●

Glossary
 C: with correction collar
 Oil: oil immersion
 WI: water immersion
 W: water dipping
 Mi: multi immersion
 Glyc: glycerin immersion
 Sil: silicone oil immersion

Cover glass thickness
 — : can be used without cover glass
 0: use without cover glass

Darkfield
 Possible with the following
 △ : universal condenser (dry) and darkfield ring
 ○ : above and darkfield condenser (dry)
 ● : darkfield condenser (oil)

Phase rings
 PHL, PH1, PH2, PH3: condenser cassette modules.
 EXT PH3, EXT PH4: external phase contrast modules for Ti2-E.

Fluorescence (UV)
 △ : possible with visible light that has a longer wavelength than the excitation light used for DAPI
 ○ : suitable
 ⊙ : recommended for best results
 340: high transmittance with an ultraviolet wavelength range of down to 340 nm

Brightfield/DIC/Polarizing/Fluorescence (visible light)
 △ : possible but not recommended
 ○ : suitable
 ⊙ : recommended for best results

Ti2-E PFS
 ● : compatible with PFS

CFI Plan Apochromat VC series

With their high degree of chromatic aberration correction and exceptional resolving power, these objectives are the perfect choice for multi-color fluorescence imaging as well as brightfield and DIC applications. In particular, axial chromatic aberration has been corrected in the short wavelength range, making these objectives highly effective for confocal applications.



CFI Plan Apochromat VC 100X Oil, 60XC WI and 20X

CFI SR series/CFI HP series/CFI SR HP series

Strict adjustment and inspection minimize the SR objective's asymmetric aberration, optimizing it for the N-SIM S super-resolution microscope. The HP objective has improved axial chromatic aberration correction and is compatible with the high-power lasers required for the fast blinks of fluorophores under the N-STORM super-resolution microscope. The SR HP objective offers both excellent optical performance and high durability against high power laser excitation, eliminating the need to switch objectives between the N-SIM S and N-STORM systems. AC objectives support the Auto Correction function of the Ti2-E inverted microscope.



CFI SR Plan Apochromat IR 60XC WI, CFI HP Plan Apochromat VC 100X Oil, CFI SR HP Apochromat TIRF 100XC Oil and CFI SR HP Plan Apochromat Lambda S 100XC Sil

Use	Model	NA	W.D. (mm)	Cover glass thickness	Brightfield	Darkfield	DIC	Phase contrast	Polarizing	Fluorescence			Ti2-E PFS
										Visible light	UV	NR	
Brightfield (CFI Plan Apo)	VC 20X	0.75	1.00	0.17	◎	○●	○		○	◎	○		●
	VC 60XC WI	1.20	0.31-0.28	0.15-0.18	◎		○	EXT PH3-60X	○	◎	◎		●
	VC 100X Oil	1.40	0.13	0.17	◎		○	EXT PH3-100X	○	◎	△		●
Super-resolution (CFI SR Plan Apo)	IR 60XC WI	1.27	0.18-0.16	0.15-0.19	◎		○	EXT PH3-60X	○	◎	○	◎	●
	IR 60XC WI	1.27	0.18-0.16	0.15-0.19	◎		○	EXT PH3-60X	○	◎	○	◎	●
Super-resolution (CFI HP Plan Apo)	VC 100X Oil	1.40	0.13	0.17	◎		○	EXT PH3-100X	○	◎	△		●
Super-resolution (CFI SR HP Plan Apo)	Lambda S 100XC Sil	1.35	0.31-0.29 (23°C) 0.30-0.28 (37°C)	0.15-0.19	◎		○		○	◎	◎		●
Super-resolution (CFI SR HP Apo)	TIRF 100XC Oil	1.49	0.16-0.10 (23°C) 0.15-0.09 (37°C)	0.13-0.19 (23°C) 0.14-0.20 (37°C)	◎		○	EXT PH4-100X	○	◎	△		●
	TIRF 100XC Oil	1.49	0.16-0.10 (23°C) 0.15-0.09 (37°C)	0.13-0.19 (23°C) 0.14-0.20 (37°C)	◎		○	EXT PH4-100X	○	◎	△		●

CFI Apochromat Lambda S series

These high-numerical-aperture (NA) objectives provide chromatic aberration correction over wavelengths ranging from violet to near-infrared and are ideal for multicolor confocal imaging. The LWD Lambda S 20XC WI lens has a chromatic aberration correction range up to infrared. The transmittance of these lenses is enhanced over a wide wavelength range by utilizing Nano Crystal Coat technology.



CFI Apochromat Lambda S 40XC WI/60X Oil and CFI Apochromat LWD Lambda S 20XC WI/40XC WI

CFI Apochromat TIRF series

These objectives boast an unprecedented NA of 1.49 (using a standard coverslip and immersion oil), the highest resolution among Nikon objectives. Correction collars enable optimization of point spread functions for varying imaging temperatures, ensuring the highest quality confocal images whether you are imaging at room temperature or 37 degrees Celsius.



CFI Apochromat TIRF 60XC Oil and 100XC Oil

Multiphoton confocal objectives

These objectives are optimized for deep tissue imaging using a multiphoton confocal microscope, thanks to their ability to correct chromatic aberration up to near-IR range. They provide both long WD and high NA, as well as high transmittance. The Glyc models have a correction collar for the refraction index of immersion liquids, and are compatible with various tissue clearing reagents that are used in neuroscience research.



CFI Plan Apochromat 10XC Glyc, CFI90 20XC Glyc and CFI75 Apochromat 25XC W 1300

Use	Model	NA	W.D. (mm)	Cover glass thickness	Brightfield	Darkfield	DIC	Phase contrast	Polarizing	Fluorescence			Ti2-E PFS
										Visible light	UV	NIR	
Confocal (CFI Apo)	LWD Lambda S 20XC WI	0.95	0.99-0.90	0.11-0.23	○	●	○		○	○	○	●	
	Lambda S 40XC WI	1.25	0.20-0.16	0.15-0.19	◎		○	EXT PH3-40X	○	◎	◎	●	
	LWD Lambda S 40XC WI	1.15	0.61-0.59	0.15-0.19	◎	●	○	EXT PH3-40X	○	◎	◎	●	
	Lambda S 60X Oil	1.40	0.14	0.17	◎		○	EXT PH3-60X	○	◎	◎	●	
Evanescent (CFI Apo)	TIRF 60XC Oil	1.49	0.16-0.10 (23°C) 0.13-0.07 (37°C)	0.13-0.19 (23°C) 0.15-0.21 (37°C)	◎		○	EXT PH4-60X	○	◎	△	●	
	TIRF 100XC Oil	1.49	0.16-0.10 (23°C) 0.15-0.09 (37°C)	0.13-0.19 (23°C) 0.14-0.20 (37°C)	◎		○	EXT PH4-100X	○	◎	△	●	
Multiphoton confocal (CFI Plan Apo)	10XC Glyc	0.50	Upright: 5.50 Inverted: 2.00	0-0.17	◎	○●				◎		◎	
Multiphoton confocal (CFI90)	20XC Glyc*	1.00	8.20	0	△**							◎	
Multiphoton confocal (CFI75 Apo)	25XC W*	1.10	2.00	0	◎	●			○	◎	○	○***	
	25XC W 1300*	1.10	2.00	0	◎	●	○		○	◎	○	○***	

* Dedicated for FN1 and Ni-E focusing nosepiece type

** Can only be used as a finder (chromatic aberration is corrected above 588 nm)

*** Near-infrared DIC

CFI Plan Fluor series

Featuring an extra-high transmission rate, especially in the ultraviolet wavelength, and flatness of field, this series is designed for fluorescence observation and imaging. These objectives can function as multi-purpose objectives for brightfield, fluorescence, simple polarizing, and DIC observations.



CFI Plan Fluor 4X, 10X, 20X, 40X, 60XC and 100X Oil



CFI Plan Fluor DLL 10X, 20X, 40X and 100X Oil

CFI Super Fluor series

This CFI Super Fluor series ensures a high transmission rate of ultraviolet wavelengths down to 340nm for fluorochromes like indo-1 and fura-2. Also, these objectives have improved S/N ratios for short wavelengths and have high NA, making the fluorescence images they produce significantly sharper and brighter.



CFI Super Fluor 10X, 20X, 40XC and 40X Oil

Use	Model	NA	W.D. (mm)	Cover glass thickness	Brightfield	Darkfield	DIC	Phase contrast	Polarizing	Fluorescence		Ti2-E PFS
										Visible light	UV	
Brightfield (CFI Plan Fluor)	4X	0.13	17.20	—	⊙				△	⊙	⊙	
	10X	0.30	16.00	0.17	⊙	△	○		○	⊙	⊙	●
	20X	0.50	2.10	0.17	⊙	●	○		○	⊙	⊙	
	20XC MI	0.75	0.51-0.35 (oil) 0.51-0.34 (glycerin) 0.49-0.33 (water)	0-0.17	⊙	●	○		○	⊙	⊙	
	40X	0.75	0.66	0.17	⊙	●	○		○	⊙	⊙	●
	40X Oil	1.30	0.24	0.17	⊙		○	EXT PH3-40X	○	⊙	⊙	●
	60XC	0.85	0.40-0.31	0.11-0.23	⊙	●	○		○	⊙	⊙	
	60XS Oil	0.50-1.25	0.22	0.17	⊙	●	○		○	⊙	⊙	
	100X Oil	1.30	0.16	0.17	⊙		○		○	⊙	⊙	●
100XS Oil	0.50-1.30	0.16	0.17	⊙	●	○		○	⊙	⊙		
Phase contrast (CFI Plan Fluor)	DL 4XF	0.13	16.50	1.20	○			⊙ PHL		○	○	
	DLL 10X	0.30	16.00	0.17	○	△		⊙ PH1		○	○	●
	DL 10XF	0.30	15.20	1.20	○	△		⊙ PH1		○	○	●
	DLL 20X	0.50	2.10	0.17	○	●		⊙ PH1		○	○	●
	DLL 40X	0.75	0.66	0.17	○	●		⊙ PH2		○	○	●
	DLL 100X Oil	1.30	0.16	0.17	○			⊙ PH3		○	○	●
	DM 40X	0.75	0.66	0.17	○	●		⊙ PH2		○	○	
BM 40X	0.75	0.66	0.17	○	●		⊙ PH2		○	○		
Apodized phase contrast (CFI Plan Fluor)	ADH 100X Oil	1.30	0.16	0.17	○			⊙ PH3		○	○	●
Brightfield (CFI Super Fluor)	4X	0.20	15.50	—	⊙				△	⊙	⊙ 340	●
	10X	0.50	1.20	0.17	⊙	●	○		△	⊙	⊙ 340	●
	20X	0.75	1.00	0.17	⊙	●	○		△	⊙	⊙ 340	●
	40XC	0.90	0.34-0.26	0.11-0.23	⊙	●	○		△	⊙	⊙ 340	
	40X Oil	1.30	0.22	0.17	⊙		○	EXT PH3-40X	△	⊙	⊙ 340	●
	100XS Oil	0.50-1.30	0.20	0.17	⊙	●			△	⊙	⊙ 340	

CFI S Plan Fluor series

The broadband multilayer coating realizes high transmittance from near-ultraviolet to near-infrared wavelengths, with superior chromatic correction. The correction collar allows these objectives to be used with a diverse range of culture vessels and specimen thicknesses. High-quality images with no aberrations can be obtained under a broad range of illumination techniques.



CFI S Plan Fluor LWD 20XC and LWD ADM 20XC



CFI S Plan Fluor ELWD 20XC, 40XC and 60XC



CFI S Plan Fluor ELWD ADM 20XC, 40XC and ADL 60XC

Nikon Advanced Modulation Contrast (NAMC) series

Nikon has developed dedicated objectives for advanced modulation contrast. Colorless and transparent samples can be observed in high relief with a plastic dish, which is not possible in DIC observation. The direction of contrast can be matched to S Plan Fluor ELWD NAMC objectives, thereby allowing optimal contrast selection for techniques like microinjection and ICSI.



CFI S Plan Fluor ELWD NAMC 20XC and 40XC



CFI NAMC 10XF, CFI LWD NAMC 20XF and 40XC

Use	Model	NA	W.D. (mm)	Cover glass thickness	Brightfield	Darkfield	DIC	Phase contrast	Polarizing	Fluorescence		Ti2-E PFS
										Visible light	UV	
Brightfield (CFI S Plan Fluor)	LWD 20XC	0.70	1.30-2.30	0-1.8	⊙	○●	○		○	⊙	⊙	●
	ELWD 20XC	0.45	8.20-6.90	0-2.00	⊙	○●	○		○	⊙	⊙	●
	ELWD 40XC	0.60	3.60-2.80	0-2.00	⊙	○●	○		○	⊙	⊙	●
	ELWD 60XC	0.70	2.60-1.80	0.10-1.30	⊙	○●	○		○	⊙	⊙	●
Apodized phase contrast (CFI S Plan Fluor)	LWD ADM 20XC	0.70	1.30-2.30	0-1.8	○	○●		⊙ PH2		○	○	●
	ELWD ADM 20XC	0.45	8.20-6.90	0-2.00	○	○●		⊙ PH1		○	○	●
	ELWD ADM 40XC	0.60	3.60-2.80	0-2.00	○	○●		⊙ PH2		○	○	●
	ELWD ADL 60XC	0.70	2.60-1.80	0.10-1.30	○	○●		⊙ PH2		○	○	●
Apodized phase contrast (CFI S Plan Fluor)	ELWD NAMC 20XC	0.45	8.20-6.90	0-2.00	○					○	○	
	ELWD NAMC 40XC	0.60	3.60-2.80	0-2.00	○					○	○	
Advanced modulation contrast (CFI)	NAMC 10XF	0.25	6.20	1.20	○					△		
	LWD NAMC 20XF	0.40	3.10	1.20	○					△		
	LWD NAMC 40XC	0.55	2.70-1.70	0-2.00	○					△		

CFI Plan Achromat series

CFI Plan Achromat series provides incredible image flatness, with chromatic aberration corrected throughout the entire visible spectrum. These objectives are suitable not only for observation but also for capturing images.



CFI Plan 1X, 2X, 4X, 10X, 20X, 40X and 100X Oil



CFI Plan DL 10X, 20X, 40X and 100X Oil

Water Dipping Series

Sharper tips and broad approach angles provide improved accessibility for manipulator control. Aberrations are corrected even in the infrared range for high-magnification objectives, making them suitable for multi-photon imaging using infrared light.



CFI Plan Fluor 10XW, CFI75 LWD 16XW, CFI Apochromat NIR 40XW/60XW and CFI Plan 100XC W

Use	Model	NA	W.D. (mm)	Cover glass thickness	Brightfield	Darkfield	DIC	Phase contrast	Polarizing	Fluorescence		Near- infrared DIC
										Visible light	UV	
Brightfield (CFI Plan)	1X	0.04	3.20	—	⊙				△	△		
	2X	0.06	7.50	—	⊙				△	△		
	4X	0.10	30.00	—	⊙				△	○		
	10X	0.25	10.50	—	⊙	△			△	○		
	20X	0.40	1.20	0.17	⊙	●			△	○		
	40X	0.65	0.56	0.17	⊙	●			△	○		
	50X Oil	0.90	0.35/0.18	-0.17	⊙	●			△	○		
	100X Oil	1.25	0.20	0.17	⊙				△	○		
	NCG 40X	0.65	0.48	0	⊙	●			△	○		
NCG 100X	0.90	1.00	0	⊙	●			△	○			
Phase contrast (CFI Plan)	DL 10X	0.25	10.50	—	○	△		⊙ PH1	△	△		
	DL 20X	0.40	1.20	0.17	○	●		⊙ PH1	△	△		
	DL 40X	0.65	0.56	0.17	○	●		⊙ PH2	△	△		
	DL 100X Oil	1.25	0.20	0.17	○			⊙ PH3	△	△		
DIC (CFI Plan Fluor)	10X W	0.30	3.50	0	⊙	△	○		○	⊙	○	
IR-DIC (CFI Apo)	NIR 40X W	0.80	3.50	0	⊙	●	○		○	⊙	△	
	NIR 60X W	1.00	2.80	0	⊙	●	○		○	⊙	○	
DIC (CFI Plan)	100XC W	1.10	2.50	0	⊙	●	○		○	⊙	○	
DIC (CFI75)	LWD 16X W *	0.80	3.00	0	⊙	●	○		○	⊙	○	

* Dedicated for FN1 and Ni-E focusing nosepiece type

CFI Achromat series

This series of objectives provide dramatic correction for chromatic aberration, spherical aberration and coma. Image flatness is significantly improved.

Nikon specifically developed ADL series for phase contrast observations by using its proprietary apodization process to improve the objective's phase ring. Imaging cell division is often impeded by unwanted halos due to the sample thickness. ADL mitigates this effect for clearer observation of mitotic events.



CFI 4X, 10X, CFI LWD 20X, CFI 40X and 100X Oil



CFI DL 10X, CFI LWD DL 20X, CFI DL 40X and 100X Oil



CFI ADL 10XF, CFI LWD ADL 20XF, 40XF and 40XC

Use	Model	NA	W.D. (mm)	Cover glass thickness	Brightfield	Darkfield	DIC	Phase contrast	Polarizing	Fluorescence		Ti2-E PFS
										Visible light	UV	
Brightfield (CFI)	4X	0.10	30.00	—	⊙				△	○		
	10X	0.25	7.00	—	⊙	△			△	○		
	LWD 20X	0.40	3.90	0.17	⊙	○●			△	○		
	40X	0.65	0.65	0.17	⊙	○●			△	○		
	LWD 40XC	0.55	2.70-1.70	0-2.00	⊙	○●			△	○		
	60X	0.80	0.30	0.17	⊙	●			△	○		
	100X Oil	1.25	0.23	0.17	⊙				△	○		
Phase contrast (CFI)	100XS Oil	0.50-1.25	0.23	0.17	⊙	○●			△	○		
	DL 10X	0.25	7.00	—	○	△		⊙ PH1	△	△		
	LWD DL 20X	0.40	3.90	0.17	○	○●		⊙ PH1	△	△		
	LWD DL 20XF	0.40	3.10	1.20	○			⊙ PH1	△	△		
	DL 40X	0.65	0.65	0.17	○	○●		⊙ PH2	△	△		
	LWD DL 40XC	0.55	2.70-1.70	0-2.00	○	○●		⊙ PH2	△	△		
	DL 100X Oil	1.25	0.23	0.17	○			⊙ PH3	△	△		
Apodized phase contrast (CFI)	BM 10X	0.25	7.00	0.70	○			⊙ PH1	△	△		
	ADL 10XF	0.25	6.20	1.20	○			⊙ PH1	△	△		
	LWD ADL 20XF	0.40	3.10	1.20	○			⊙ PH1	△	△		
	LWD ADL 40XF	0.55	2.10	1.20	○			⊙ PH1	△	△		
	LWD ADL 40XC	0.55	2.70-1.70	0-2.00	○	○●		⊙ PH2	△	△		

Specifications and equipment are subject to change without any notice or obligation on the part of the manufacturer. December 2018
© 2018 NIKON CORPORATION



WARNING

TO ENSURE CORRECT USAGE, READ THE CORRESPONDING MANUALS CAREFULLY BEFORE USING YOUR EQUIPMENT.

Company names and product names appearing in this brochure are their registered trademarks or trademarks.
N.B. Export of the products* in this brochure is controlled under the Japanese Foreign Exchange and Foreign Trade Law. Appropriate export procedure shall be required in case of export from Japan.
*Products: Hardware and its technical information (including software)



NIKON CORPORATION

Shinagawa Intercity Tower C, 2-15-3, Konan, Minato-ku, Tokyo 108-6290, Japan
phone: +81-3-6433-3705 fax: +81-3-6433-3785
<https://www.healthcare.nikon.com/>

ISO 14001 Certified
for NIKON CORPORATION

NIKON INSTRUMENTS INC.

1300 Walt Whitman Road, Melville, N.Y. 11747-3064, U.S.A.
phone: +1-631-547-8500; +1-800-52-NIKON (within the U.S.A. only)
fax: +1-631-547-0306

<https://www.microscope.healthcare.nikon.com/>

NIKON INSTRUMENTS EUROPE B.V.

Tripolis 100, Burgervveeshuispad 101, 1076 ER Amsterdam, The Netherlands
phone: +31-20-7099-000 fax: +31-20-7099-298

https://www.microscope.healthcare.nikon.com/en_EU/

NIKON INSTRUMENTS (SHANGHAI) CO., LTD.

CHINA phone: +86-21-6841-2050 fax: +86-21-6841-2060
(Beijing branch) phone: +86-10-5831-2028 fax: +86-10-5831-2026
(Guangzhou branch) phone: +86-20-3882-0550 fax: +86-20-3882-0580

NIKON CANADA INC.

CANADA phone: +1-905-602-9676 fax: +1-905-602-9953

NIKON FRANCE S.A.S.

FRANCE phone: +33-1-4516-45-16 fax: +33-1-4516-45-55

NIKON GMBH

GERMANY phone: +49-211-941-42-20 fax: +49-211-941-43-22

NIKON INSTRUMENTS S.p.A.

ITALY phone: +39-55-300-96-01 fax: +39-55-30-09-93

NIKON GMBH SWITZERLAND

SWITZERLAND phone: +41-43-277-28-67 fax: +41-43-277-28-61

NIKON UK LTD.

UNITED KINGDOM phone: +44-208-247-1717 fax: +44-208-541-4584

NIKON CEE GMBH

AUSTRIA phone: +43-1-972-6111 fax: +43-1-972-611-140

NIKON SINGAPORE PTE LTD

SINGAPORE phone: +65-6559-3651 fax: +65-6559-3668

NIKON INSTRUMENTS KOREA CO., LTD.

KOREA phone: +82-2-2186-8400 fax: +82-2-555-4415

Официальный дистрибьютор Nikon
в России и странах СНГ



Россия, 127055, г. Москва, ул. Новолесная, д. 2

тел.: +7 (495) 223-40-00 факс: +7 (495) 223-40-01

<http://www.tokyo-boeki.ru> email: systems@tokyo-boeki.ru